

Table Top Small Refrigerated Centrifuge Temperature Controlled Custom

Our Product Introduction

for more products please visit us on mkecentrifuge.com

Basic Information

- Place of Origin: CHINA
- Brand Name: MKE
- Certification: CE,ISO13485,ISO9001,SFDA
- Model Number: MI6R (old model:KCL-60R)
- Minimum Order Quantity: 1
- Price: Negotiable
- Packaging Details: Wooden Box
- Delivery Time: 7-15 work days
- Payment Terms: L/C,D/A,D/P,T/T,Western Union,MoneyGram
- Supply Ability: 100 set/month




Product Specification

- Applications: Blood Collection Trolleys, Blood Collection Houses
- Max Speed: 6000rpm
- Max RCF: 4860xg
- Max Capacity: 4*100ml Swing Bucket Rotor
- Matched Rotors: Fixed Angle Rotors, Swing Out Rotor
- Display: LCD, Touch Five-inch Colour Screen Display
- Warranty: 1 Year, 3 Years, 5 Years, Free Within Warranty, Free Within Warranty
- Highlight: Tabletop small refrigerated centrifuge, Custom small refrigerated centrifuge, Custom table top refrigerated centrifuge



More Images

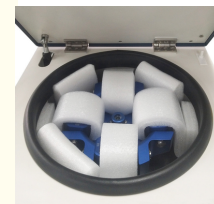

Product Description

Blood collection trolleys or houses may be used by a variety of personnel, so microcentrifuge KCL-60R designed for this purpose are typically fast and efficient, with high spin speeds and short cycle times, easy to use and intuitive. It has pre-programmed settings that make it easy to spin the blood correctly. It process blood quickly to ensure that it remains viable for medical use. Equipped with 4 x 100ml horizontal rotor, optional nucleic acid centrifuge tubes, vacuum blood collection tubes, 15ml, 10ml, 8ml, 5ml, 2ml and other centrifuge tube adapters.

Design

- Metal housing, metal lid, and stainless steel centrifuge chamber are autoclave-able and resistant detergents and disinfection liquids.
- Observations window on lid for speed calibration and noise measurement.
- Microprocessor control with frequency conversion A.C. motor drive.
- Adopt high efficiency environmental protection refrigeration system, no environmental pollution, high precision temperature control.
- LCD touch screen with simple&clear interface, display all parameters of the instrument operation, the

Our Product Introduction

user can set the centrifugal speed and centrifugal time arbitrarily.

- Centrifuge operating system with a password to prevent misuse, external USB interface, can track and trace more than 1000 sets of centrifugal data.
- Light weight, small size, simple, beautiful and generous appearance, unique suction cup feet can firmly grasp the placement surface, even when driving without worrying about centrifuge movement.
- Four suction cup feet can be adjusted high and low, while the centrifuge built-in level, blood collection vehicle in the case of uneven ground can still adjust the centrifuge to a stable state of use, the best centrifugal effect, the use of safer.
- Equipped with 4 x 100ml horizontal rotor, optional nucleic acid centrifuge tubes, vacuum blood collection tubes, 15ml, 10ml, 8ml, 5ml, 2ml and other centrifuge tube adapters.



User friendly

- CFC free refrigeration system, precooling rotors function during standby, precise temperature control up to -9°C .
- Programmable and storable, stores up to 1000 programs and provides direct user access with one keystroke save. The stored parameter: combination remains in the memory even after the centrifugation has been switched off.
- Automatic lid-locking and holding during rotor run
- Easily exchangeable rotors
- Automatic troubleshooting with indication of error code and alarm on the screen.
- "Mute" patent technology. With low noise "patent" certificate, the highest speed operation noise is less than 50 decibels.



Safety

- Produced according to international safety regulations (e.g. IEC 61010)
- Safety Electrical lid interlock, door can not open when rotor is running and machine won't run when lid open and error code appears.
- Emergency lid-lock release when power failure.
- Gas hinge for lid drop protection.
- Safety protections are made for the over speed, over heating (motor, chamber), imbalance etc to ensure the machine safety.



WHY CHOOSE US

- **Direct from the Manufacturer:** Enjoy competitive pricing on all our products.
- **Certified Excellence:** Experience premium quality with our CE and ISO certified products.
- **Customization at Its Best:** Benefit from our robust OEM and ODM capabilities tailored to your specific needs.
- **Reliable Production and Delivery:** Trust us for guaranteed on-time production and delivery of your orders.
- **Comprehensive Support:** Rely on our comprehensive after-sales service, including complimentary product accessories.
- **Expertise You Can Trust:** We stand behind our 15 years of experience in R&D and production of medical centrifuge machines.



Our Factory Introduction

We have a standardized production workshop operated by a dedicated team that ensures quality in every laboratory centrifuge and lab equipment we produce. Our R&D division includes design engineers, technology specialists, qualified experts, and process control managers. Following a strict quality control system, our specialists oversee the entire production process, from design and R&D to testing, production, quality inspection, and packaging.

- Specialize in centrifuges.
- With decades of experience, we've exported to 80+ countries.
- Proven Reputation: SGS and Intertek inspected, underlining our reliability.
- Manufacturing Capabilities: 6 floors, 3000 sqm each, dedicated to manufacturing and quality control.




MKE Exhibition




Packaging & Delivery

MKE centrifuges are securely packaged in wooden cases or carton boxes, selected according to the weight and size of the goods. This ensures they are delivered in excellent condition to customers worldwide.

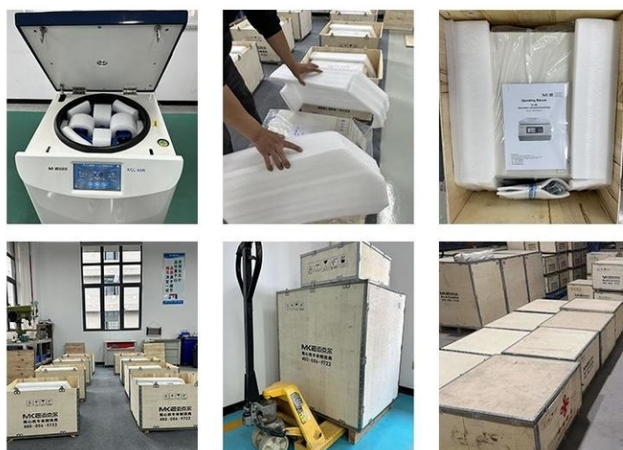

After-Sale WARRANTY

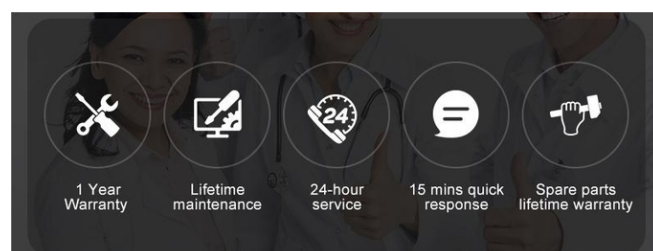

FAQ

Q: How long is the production lead time?

A: Within 7-15 workdays after receiving deposit.

Q: What's your warranty?

A: All MKE products have ONE year warranty, five-year motor warranty. We won't charge clients any money for repairing the machine during one-year warranty time, the freight is also on us.

Q: How about the sample of MKE Centrifuge?

A: We could send you the samples but the freight and the sample are charged.

Q: Which kinds of Certificate would you have?

A: ISO9001, ISO13485, CE

Q: Do you offer centrifuges with voltage and plug options suitable for my country/region?

A: Yes, we provide centrifuges with various voltage and plug options to meet different country/region requirements. Please let us know your specific needs, and we will ensure compatibility.

Q: Are there customized options available for specific applications?

A: Absolutely! We offer customization options to meet specific application needs. Let us know your requirements, and we will tailor the centrifuge accordingly.



Professional
Centrifuge
Manufacturer

Hunan Michael Laboratory Instrument Co., Ltd.

+86-17775883371

Shirley@mke-lab.com

mkecentrifuge.com

MKE Building, Sanhe Zhihui Industrial Park, Guyuan Road No. 426, Changsha High-tech Development Zone, Hunan, China