

Benchtop Clinical Centrifuge With Swing Rotor Efficiently Process 12 Cell Slides

Our Product Introduction

Basic Information

- Place of Origin: CHINA
- Brand Name: MKE
- Certification: CE,ISO13485,ISO9001,SFDA
- Model Number: XJ-6
- Minimum Order Quantity: 1
- Price: Negotiable
- Packaging Details: Wooden Box
- Delivery Time: 7-15 work days
- Payment Terms: L/C,D/A,D/P,T/T,Western Union,MoneyGram
- Supply Ability: 100 set/month



Product Specification

- Product Type: Benchtop Cell Cytospin Centrifuge
- Applications: Ideal For The Centrifugation Of Samples In Clinical, Chemistry, Cytology, Research And Industrial Laboratories.
- MaxSpeed: 1200rpm
- MaxRCF: 241xg
- MaxCapacity: 12 Cell Slide Clips
- Noise Level: ≤55dB(A)
- MatchedRotors: Swing Rotor
- Display: LED
- Warranty: 1 Year Warranty For The Whole Machine, 3 Years Warranty For The Compressor And 5 Years Warranty For The Motor. Free Replacement Parts And Shipping Within Warranty.
- Highlight: **Benchtop Clinical Centrifuge, 12 Cell Slides Clinical Centrifuge, Swing Rotor Clinical Centrifuge**



More Images



for more products please visit us on mkecentrifuge.com

Product Description

Description

The XJ-6 Cell Smear Centrifuge then processes diagnostic components like cells and bacteria into a flat, uniform layer, free from overlaps.

These ultra-thin layers adhere to glass slides, ready for microscopic examination or computer analysis.

With the capacity to process 12 samples in just 5 minutes, it combines user-friendly operation, eco-friendly design, and consistent cell layering, even on a clear background.

Samples collected via cervical brushing are purified using gradient centrifugation, effectively removing impurities.



Specifications:

Speed range(rpm)	0-1200rpm adjustable	Noise(dBA)	≤55 dB(A)
Control	Microprocessor control	Motor	Inverter motor direct drive
Max Capacity(ml)	12 cell slide clips	Height of Cell Slide Clips	50mm
Time Range	0-99min	Dimension(mm)	440*400*230mm
RCF range (xg)	0-241xg adjustable	Net weight(without rotor)	22KG
Power supply	AC220V,50HZ,5A,60W		

Features

The centrifuge boasts a durable metal housing, metal lid, and stainless steel chamber, all of which are autoclave-able and resistant to detergents and disinfectants.

Precision is ensured with a microprocessor-controlled system powered by a frequency conversion A.C. motor.


An LED display provides clear readings of speed, timer, RCF, rotor programs, acceleration, deceleration,

and error codes for easy monitoring.

A double-layer sealing ring secures the lid and instrument, ensuring a tight and reliable seal.

Advanced bio-safety features include a gas-tight rotor head sterilizable at high temperatures and a UV sterilization function within the chamber, automatically activating during startup.

Matched Rotor

No.	Rotor	Capacity	Max. speed	Max. RCF	Note
NO.1 	Swing rotor	12 cell slide clips 2ml/clip	1200r/min	241xg	Swing rotor of settlement type



WHY CHOOSE US

- **Direct from the Manufacturer:** Enjoy competitive pricing on all our products.
- **Certified Excellence:** Experience premium quality with our CE and ISO certified products.
- **Customization at Its Best:** Benefit from our robust OEM and ODM capabilities tailored to your specific needs.
- **Reliable Production and Delivery:** Trust us for guaranteed on-time production and delivery of your orders.

- **Comprehensive Support:** Rely on our comprehensive after-sales service, including complimentary product accessories.
- **Expertise You Can Trust:** We stand behind our 15 years of experience in R&D and production of medical centrifuge machines.



Our Factory Introduction

We have a standardized production workshop operated by a dedicated team that ensures quality in every laboratory centrifuge and lab equipment we produce. Our R&D division includes design engineers, technology specialists, qualified experts, and process control managers. Following a strict quality control system, our specialists oversee the entire production process, from design and R&D to testing, production, quality inspection, and packaging.

- Specialize in centrifuges.
- With decades of experience, we've exported to 80+ countries.
- Proven Reputation: SGS and Intertek inspected, underlining our reliability.
- Manufacturing Capabilities: 6 floors, 3000 sqm each, dedicated to manufacturing and quality control.




MKE Exhibition




Packaging & Delivery

MKE centrifuges are securely packaged in wooden cases or carton boxes, selected according to the weight and size of the goods. This ensures they are delivered in excellent condition to customers worldwide.

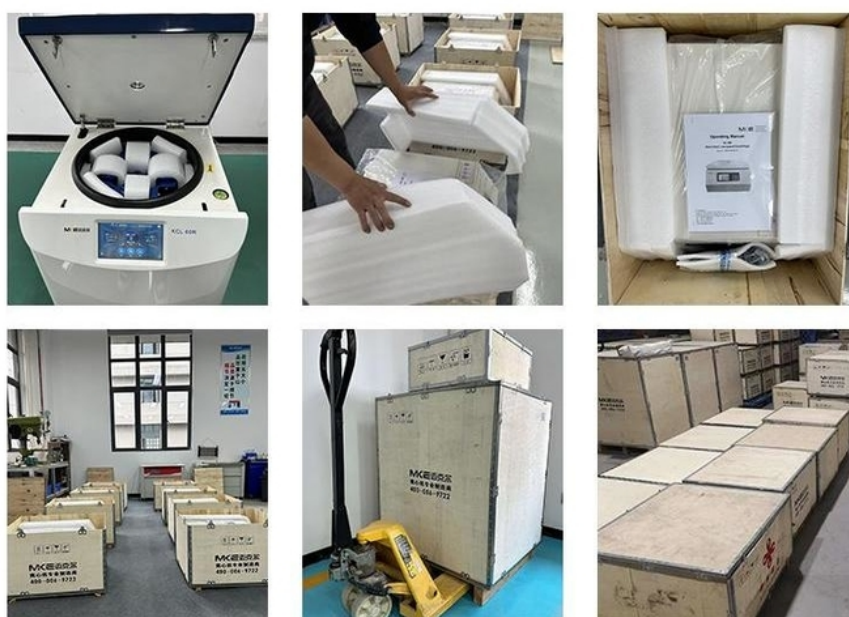

After-Sale WARRANTY

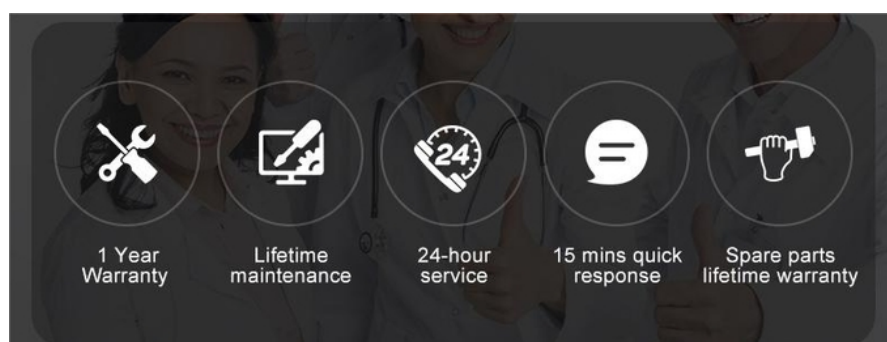

FAQ

Q: How long is the production lead time?

A: Within 7-15 workdays after receiving deposit.

Q: What's your warranty?

A: All MKE products have ONE year warranty, five-year motor warranty. We won't charge clients any money for repairing the machine during one-year warranty time, the freight is also on us.

Q: How about the sample of MKE Centrifuge?

A: We could send you the samples but the freight and the sample are charged.

Q: Which kinds of Certificate would you have?

A: ISO9001, ISO13485, CE

Q: Do you offer centrifuges with voltage and plug options suitable for my country/region?

A: Yes, we provide centrifuges with various voltage and plug options to meet different country/region requirements. Please let us know your specific needs, and we will ensure compatibility.

Q: Are there customized options available for specific applications?

A: Absolutely! We offer customization options to meet specific application needs. Let us know your requirements, and we will tailor the centrifuge accordingly.



Professional
Centrifuge
Manufacturer

Hunan Michael Laboratory Instrument Co.,Ltd.



+86-17775883371



Shirley@mke-lab.com



mkecentrifuge.com

MKE Building, Sanhe Zhihui Industrial Park, Guyuan Road No. 426, Changsha High-tech Development Zone,
Hunan, China